

# Stimulus calibration

---

## [Stimulus calibration](#)

### [Overview](#)

#### [Principle](#)

#### [Calibration in 4 steps:](#)

### [Hardware](#)

#### [Photometer](#)

#### [Arduino UNO Broad assembling](#)

#### [Assembling](#)

### [Software](#)

#### [Preparation for calibration](#)

##### [Part A: SimpleAIRecorder.ino](#)

##### [Part B. Calibration.ino Might be different in the other setup](#)

##### [Part C. iPython notebook](#)

### [Run Calibration Measurement](#)

### [FAQ](#)

[Q: The readout of photometer is not stable or too small](#)

[Q: The photometer's readout is out of range](#)

[Q: Photometer turned itself off..](#)

[Q: Takes forever to upload the Arduino sketches](#)

[Q: Arduino returns the error "access is denied...." when uploading the sketches to the port](#)

[Q: Random weird symbols instead of numbers displayed on the SimpleAIRecorder serial monitor](#)

[Q: A1 inputs shown in the serial monitor are out of range \(e.g. A1 input = 1023\) or stay at zero:](#)

[Q: Error: "No module named 'simple\\_ai\\_recorder'" when running the 2nd code block](#)

[Q: Error: "No module named 'serial'" when running the 2nd code block](#)

[Q: Error: "\[Errno 2\] No such file or directory" when running the 6/7/13th code block](#)

[Q: No messages or plots returned after starting the calibration measurement \(11th code block\)](#)

[Q: 'Error: Link not open' when running the 11th block](#)

[Q: 'Error: list index out of range' when running the 12th block, no errors returned in the 11th block](#)

[Q: 'Error: Less power levels recorded ...' in 12th block](#)

## Overview

---

### Principle

- What are we calibrating : **Power of LED emitted light** to make the **photoisomerization rate of mouse M and S opsin** to **green** and **blue light stimulus** respectively be **the same**
- How do we calibrate :
  - Measure the **input power** and the **light intensity** ( **I/O** ) of LEDs
  - Calculate **the relationships between the input power and the photo-isomerization rate** of M and S opsins for green and blue LEDs from the LED **I/O**
  - **Calibrate the input power** accordingly
- See [Calculating R star](#) and the **ipython notebook** **stimulator\_calibration\_v4.ipynb** in the [documentation folder](#) for more explanation.

## Calibration in 4 steps:

- Collect the information of the **spectrum of LED, LED filter, dichroic beam splitter**
- Assemble the hardware for LED **I/O** measurement → see [Hardware](#)
- Upload the **ipython notebook** and **Arduino sketches** for LED measurement and plotting the relationships between the input power and the photoisomerization rate → see [Software](#)
- Run the **ipython notebook** and get the plots, and adjust the input power accordingly → see [Run Calibration Measurement](#)
- Usually the calibration can be done within 5-10 mins, read [FAQ](#) if you got some unexpected problems

**Note:** I have only tried the calibration on Setup 1. The procedures **might need to be done differently for the other setups** are **highlighted**.

---

## Hardware


---

### Photometer

- Power Cable for photometer (**Label 1**):



- Others:


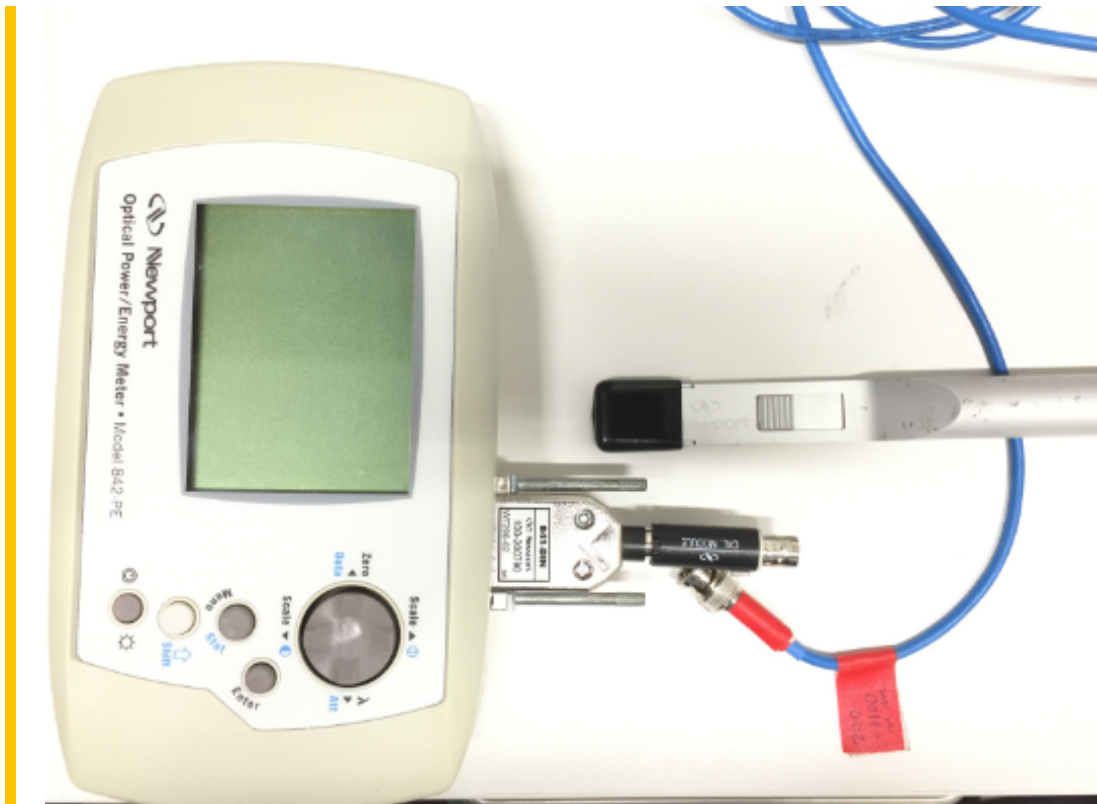


## Arduino UNO Broad assembling



- **UNO** Broad:



- Cable **[2]** (for trigger input) with two thin cables **A0** and **GND** :

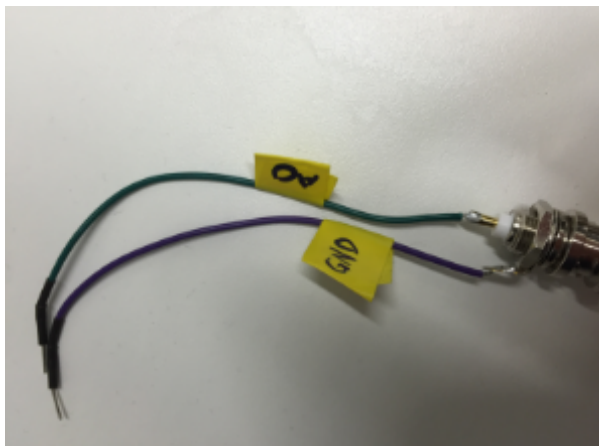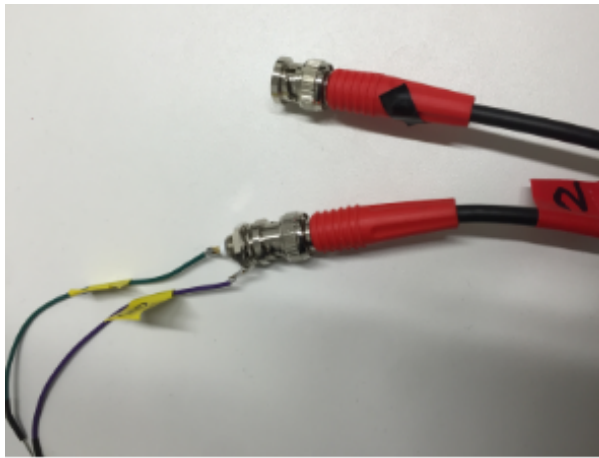

- Cable [3] for photometer input with two thin cables A1 and GND :

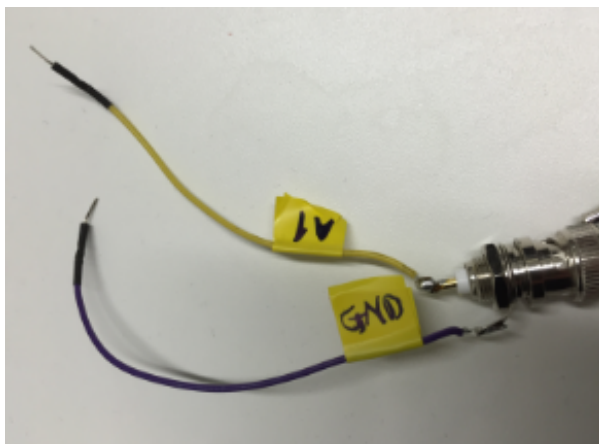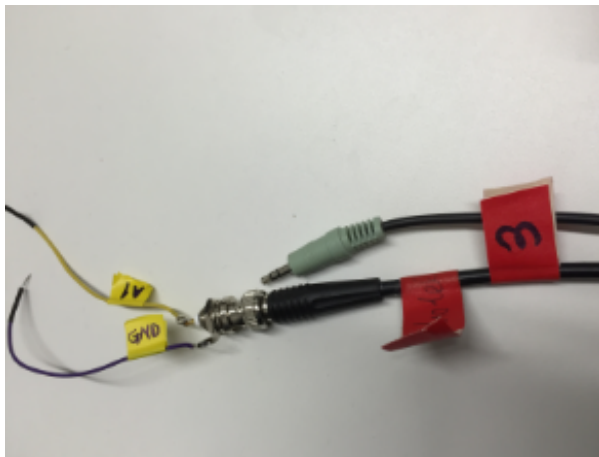


- **USB** cable:



## Assembling

- Connect **USB** to the port on Arduino board and on PC
- Connect the thin cables with **GND** label to the **GND** ports on **UNO** (GND: Ground)
- Connect the thin cables **A0** (on Cable [2]) and **A1** (on Cable [3]) to the **A0** and **A1** ports on **UNO** respectively:

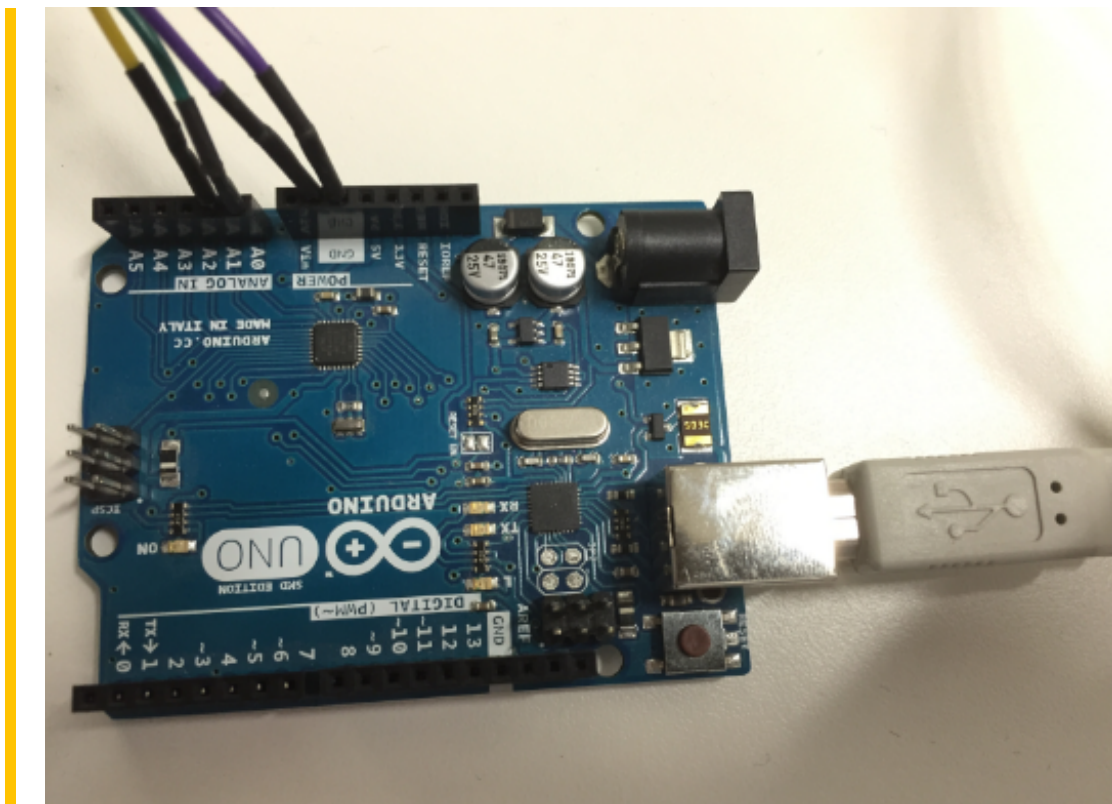


- Connect the other end of Cable **[2]** to the port for trigger channel


**I have only labeled the trigger channel port for Setup 1**, you need to find out the right port for the other setups

- Assemble the photometer: Plug in Cable **[1]** and **[3]** :




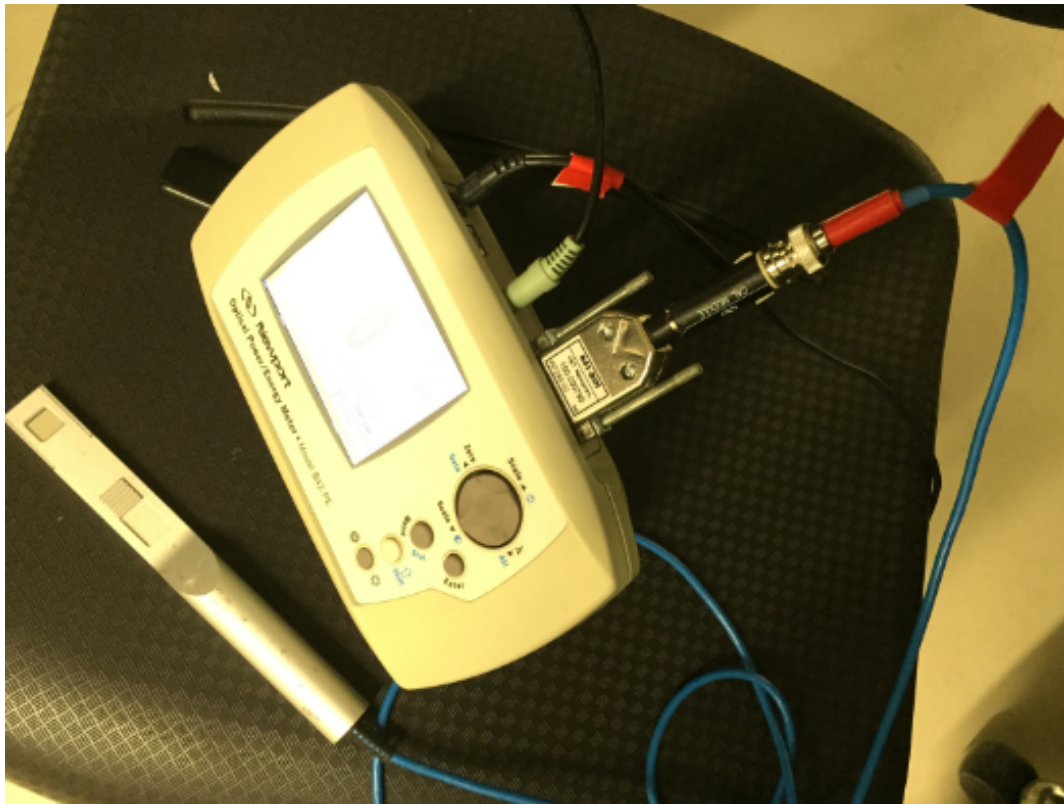

- Switch on the photometer, check the **data sampling rate** (set to 20 pts/sec, can be higher):

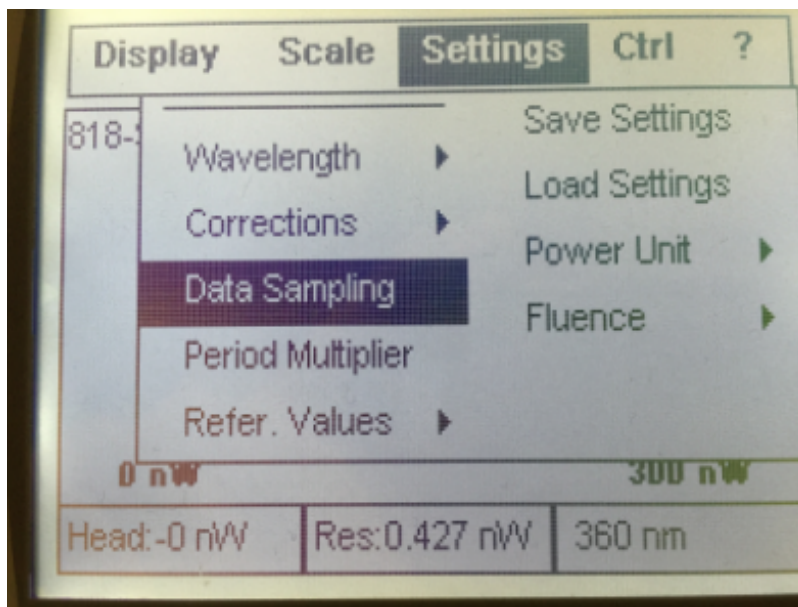

- Turn on the **analog output**:

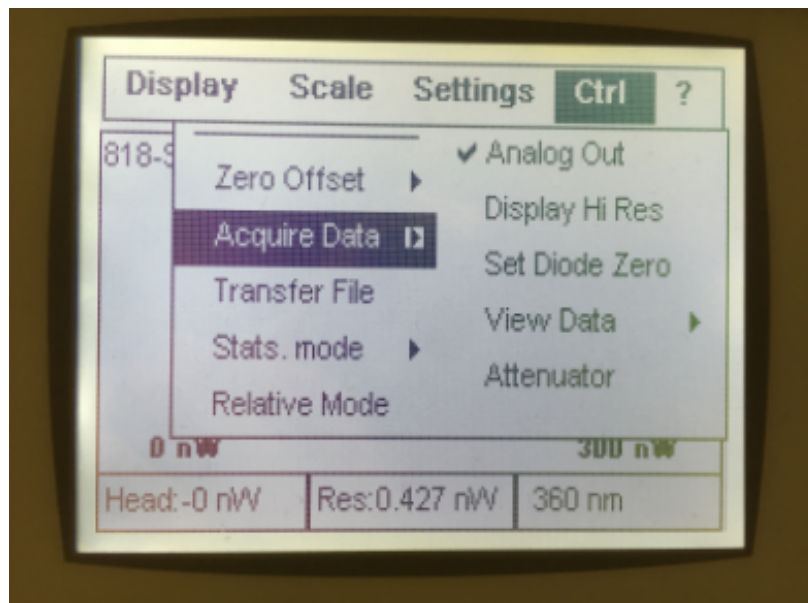


- Turn on green and blue LED, **turn off infrared LED** to avoid any interference :





Might be different in the other setup

- **Adjust the condenser** (below the microscope chamber) so the photometer can get max light input:

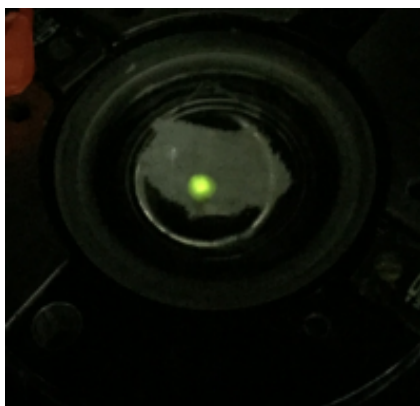

Might be different in the other setup

- **Expose the photometer sensor and place the photometer sensor area in the centre of the chamber.**  
**Measure in dark** to avoid the interference of other light sources:

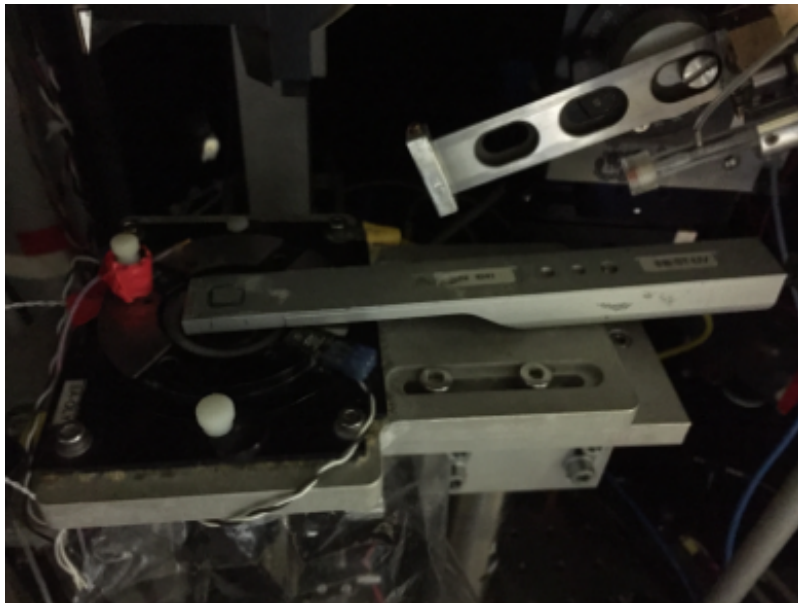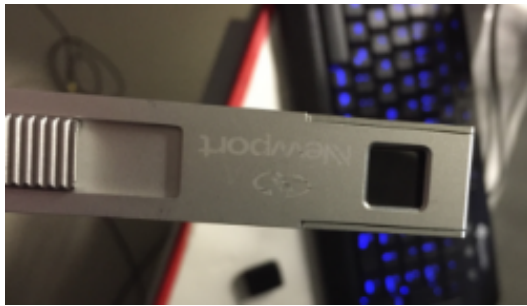

---

## Software

---

- General: **Igor**, **jupyter notebook**, **python**, **pyserial python module**
- The iPython notebook **stimulator\_calibration\_v4.ipynb** in the ["softwares"](#) folder
- **SimpleAIRRecorder.ino** in the ["arduino"](#) folder
- **Calibration.ino** in the ["calibration"](#) folder (if your LED is controlled by Arduino broad)

### Preparation for calibration

#### Part A: SimpleAIRRecorder.ino

- This sketch file returns the trigger time and light intensity data to PC so the iPython notebook can access the data for calibration.
- Open and upload the **SimpleAIRRecorder.ino** sketch to **UNO** board in Arduino IDE, **before you upload**, make sure:


Board type: **UNO**

---




Port: **COM...** (Arduino Uno) (usually **COM9 (Arduino Uno)** but can be changed

sometimes)

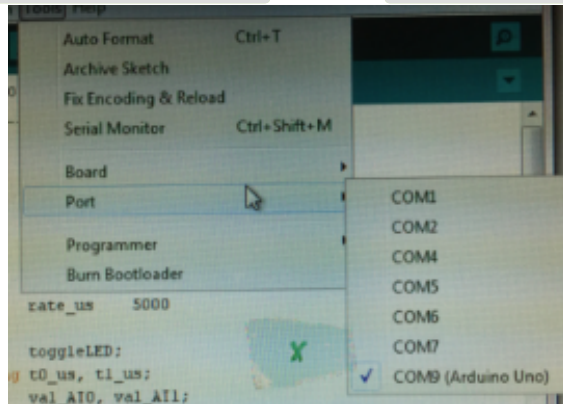

- Open the **serial monitor** (the rightmost icon on the image above) and switch baud to **115200 baud**, you should see something like this:

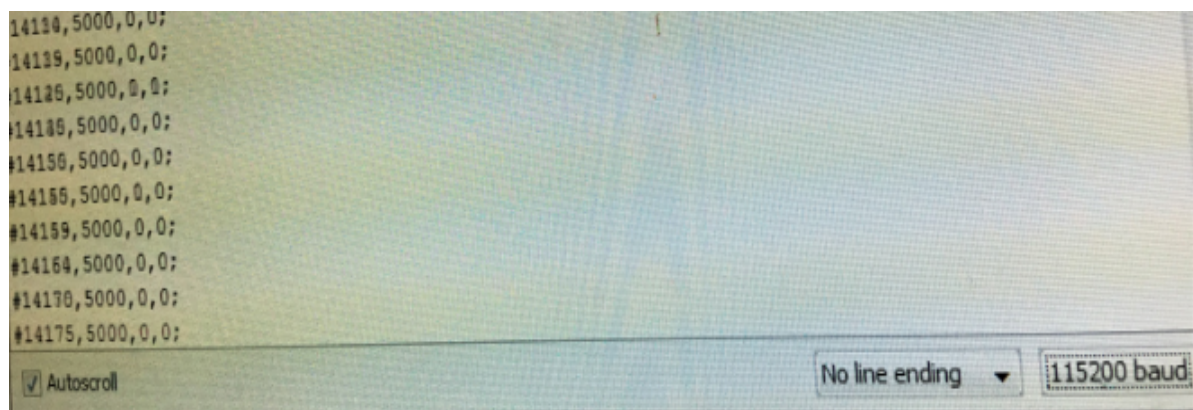

- For each line from left to right: **time (ms)**, **measuring rate ( $\mu$ s)**, **A0 (trigger) input**, **A1 (photometer) input**
- Switch on the photometer and see if the A1 input returns nonzero number, if not, check if you assembled everything correctly.

## Part B. Calibration.ino **Might be different in the other setup**

- This sketch file contains **the calibration stimuli for green and blue LEDs** which allow to measure the power for a series of light levels applied to the different LEDs in sequence
- Switch on the **microcontroller** and open **Igor**, set a configuration and start to scan (**otherwise the LEDs won't response**)
- Upload the **Calibration.ino** sketch to **LED controlling Arduino board**, again **check the board type and the port**

In **Setup 1** the board type is **Duemilanove** and the port name is **COM6**.

- Open the serial monitor for the **Calibration.ino** sketch, set **green or blue LEDs to maximum** (commands = **2** and **4** respectively) in turns to **adjust the measuring range of the photometer to a reasonable range**:
  - Make sure you change the **measuring wavelength of the photometer** ( **$\lambda$**  on the photometer panel) to the corresponding green and blue light wavelength and **'Zero'** the photometer before every measurement.

- **Make sure you set the LED you are not measuring to minimum level** ( **1** for green LED, **3** for blue)
- Play a white flash stimuli (command = **8**) and see if the change of the **A0** and **A1** values returned in the **serial monitor for SimpleAIRRecorder sketch** match the stimuli.

### Part C. iPython notebook

- Type **jupyter notebook** in **Windows Powershell**, which should open the **Home page** of the notebook. Upload the iPython notebook and open it in a new tab.
- The notebook contains **the script for calibration** and also **detailed explanation of the reasons and principles of the calibration** which I would highly recommend to read through before your first calibration.
- The following parameters are required by the script and you need to know them all before the calibration:
  - **LEDs and LED filters spectrum**

**For Setup 1**, a narrow filter is used therefore only the peak and the bandwidth of the LED/filter spectrum )

**For the other setup**, .txt files of the spectrum can be found in the [data](#) folder (e.g. "F73-063\_z400-580-890.txt")
  - **Spectrum of dichroic beam splitter:**

**For setup 1 not needed**, because the spectrum of the beam splitter is mostly overlapped with the spectrum of the two LEDs
  - **Mouse M and S Opsin spectrum:** "mouse\_cone\_opsins.txt" in the [software](#) folder
- The paths: **script\_path** and **pathData** In the **2nd and 3rd python code block** of the notebook are the path where your **ipython notebook, mouse opsin spectrum and LED/filters spectrum files** are located respectively. **pathData** is also the path to **save all the calibration measurements**
- In the **10th code block (Read calibration data)**, change the **comPortName** if it does not match with the port name of the **UNO Arduino board**.
- The **calFileName** is the **file name for the next calibration measurement**, **change it every time before you start a new measurement, otherwise no data will be saved.**

## Run Calibration Measurement

1. Run the **11th code block** in the notebook, and you should see the following:

```

Illuminated area is 1.963 mm2
Prepare to measure LED `green` and press Enter to continue...

```

2. **Set the measuring wavelength of the photometer to the peak of the corresponding LED spectrum** (e.g. in setup 1, measuring wavelength for green LED should be set to 578)
3. In the **serial monitor window for the Calibration.ino sketch**, type **g** which is the **command for green LED calibration measurement**.
4. Wait until seen this:

```

Illuminated area is 1.963 mm2
Prepare to measure LED `green` and press Enter to continue...
Opened serial port COM6 at 115200 baud
Sampling rate is 5.001 ms
56.000 s duration = 11198 samples
100% done
SUCCESS
11198 data points recorded
Rate = 5.000 +/- 0.001 ms
Closed serial port COM6
Prepare to measure LED `blue` and press Enter to continue...

```

5. Repeat **Step 2 and 3** for the other LED (type **b** for blue LED calibration measurement).
6. Run the rest of the script until you get **the plot of photoisomerization rates for each LED/filter vs. photoreceptor combination**

---

## FAQ

---

### Q: The readout of photometer is not stable or too small

- a. Make sure the cables are not broken.
- b. examine **the connection between the photometer calibration module and the photometer device** (where your photometer sensor connect to the device)
- c. Adjust the condenser and the position of the photometer sensor to make sure the sensor can receive the maximum output
- d. Adjust the scale of the photometer to a reasonable range
- e. Make sure you **take off the cap and pull down the filter to really expose the photometer sensor**
- f. Check if the hardwares are wired correctly (e.g. **A0**, **A1**, **GND** really connect to the right ports). Make sure you don't shunt anything

### Q: The photometer's readout is out of range

- a. Make sure the sensor is placed in dark so **no other light source** can interfere the photometer measurement.
- b. Adjust the scale of the photometer to a reasonable range
- c. Check the wiring and make sure you don't shunt anything
- d. Check if **the IR-LED has been turned off**

### Q: Photometer turned itself off..

- Plug in the power cable ... (one stupid mistake I repeatedly made)

### Q: Takes forever to upload the Arduino sketches

- a. Check if you select the **correct broad type** and **ports** on the Arduino IDE
- b. Check your **firewall settings** and open and upload with the administrator permission
- b. Check if you can successfully compile your sketches

**Q: Arduino returns the error "access is denied...." when uploading the sketches to the port**

The Arduino broad has been occupied by some other processes (usually the ipython notebook and the python functions it called), you must release it before uploading the sketch.

- a. **Restart your ipython notebook**
- b. **End the python processes**
- c. Reboot your Arduino broad by re-plugging in the USB cable
- d. Restart the PC

**Q: Random weird symbols instead of numbers displayed on the SimpleAIRRecorder serial monitor**

Switch baud to **115200 baud**

**Q: A1 inputs shown in the serial monitor are out of range (e.g. A1 input = 1023) or stay at zero:**

**A1 inputs should be in the range from 0 to 999.**

- a. Check if the **A1 / A0** cables are **mistakenly plugged into the wrong port** on the **UNO** Arduino broad.
- b. Check if the **analog output** of the photometer has not been turned on
- c. Make sure Cable **[2]** is connected to the photometer device.
- d. Make sure you have turned on the LEDs (**please don't directly stared at the blue and UV LEDs**)
- d. **Turn on the microcontroller, open Igor and start recording**

**Q: Error: "No module named 'simple\_ai\_recorder'" when running the 2nd code block**

Check if you can find **simple\_ai\_recorder.py** under the **script\_path** you defined

**Q: Error: "No module named 'serial'" when running the 2nd code block**

Open **Windows powershell**, type **pip install pyserial** and return, wait until the installation finished

**Q: Error: "[Errno 2] No such file or directory" when running the 6/7/13th code block**

Check if the **"mouse\_cone\_opsins.txt"** and your **spectrum files** is under the **script\_path** and the **pathdata** you defined

**Q: No messages or plots returned after starting the calibration measurement (11th code block)**

This can be due to many reasons, here I picked the most common ones:

- a. **comPortName** does not match with the port name of the **UNO** Arduino broad.
- b. The current python kernel is occupied by some other processes, **restart the kernel**.

c. Check the **A0** / **A1** / **GND** connection and make sure you don't shunt anything

**Q: 'Error: Link not open' when running the 11th block**

Usually caused by the connection problems between the **Uno** board, the photometer and the PC. Check if you assemble the hardware correctly and you don't shunt anything

**Q: 'Error: list index out of range' when running the 12th block, no errors returned in the 11th block**

Check if the **A1** / **A0** cables are **mistakenly plugged into the wrong port** on the **UNO** Arduino board.

**Q: 'Error: Less power levels recorded ...' in 12th block**

Check **all the parameters in the 10th code block**, e.g. **nLevels** if it's smaller or greater than the number of levels you measured